



# Semrock Raman Optical Filter Catalog


## HIGH-PERFORMANCE RAMAN SPECTROSCOPY SOLUTIONS

IDEX Health & Science | Semrock offers the widest selection of Raman spectroscopy-related long pass and dichroic filters available, with edge wavelengths from 224 to 1550 nm. They offer the highest transmission and the lowest ripple, and are available with a range of edge steepness values, from the highest available to lower values, at a range of prices. These filters can be paired with MaxLine® laser line filters for the best suppression of the laser background. Now, you can see the weakest signals closer to the laser line than ever before.

**DISCOVER  
WHAT'S  
INSIDE:**



Latest Raman  
Updates

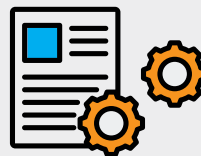

Technical  
Notes



Optical Filter  
Product  
Information




Intelligent Solutions for Life®  
Fluidics | Optics | Consumables | Assemblies

## TECHNICAL NOTE

# Filter Types for Raman Spectroscopy Applications

Raman spectroscopy is widely used today for applications ranging from industrial process control to laboratory research to bio/chemical defense measures. Industries that benefit from this highly specific analysis technique include the chemical, polymer, pharmaceutical, semiconductor, gemology, and medical fields. In Raman spectroscopy, an intense laser beam is used to create Raman (inelastic) scattered light from a sample under test. The Raman "fingerprint" is measured by a dispersive or Fourier Transform spectrometer.

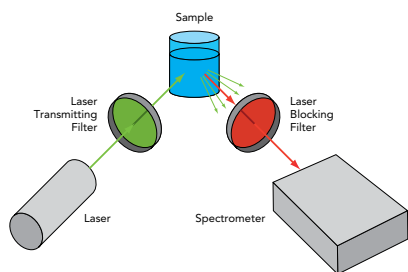


There are three basic types of Raman instrumentation. Raman microscopes, also called micro-Raman spectrophotometers, are larger-scale laboratory analytical instruments for making fast, high-accuracy Raman measurements on very small, specific sample areas. Traditional laboratory Raman spectrometers are primarily used for R&D applications, and range from "home-built" to flexible commercial systems that offer a variety of laser sources, means for holding solid and liquid samples, and different filter and spectrometer types. Finally, a rapidly emerging class of Raman instrumentation is the Raman micro-probe analyzer. These complete, compact and often portable systems are ideal for use in the field or in tight manufacturing and process environments. They utilize a remote probe tip that contains optical filters and lenses, connected to the main unit via optical fiber.

Optical filters are critical components in Raman spectroscopy systems to prevent all undesired light from reaching the spectrometer and swamping the relatively weak Raman signal. Laser Transmitting Filters inserted between the laser and the sample block all undesired light from the laser (such as broadband spontaneous emission or plasma lines) as well as any Raman scattering or fluorescence generated between the laser and the sample (as in a fiber micro-probe system). Laser Blocking Filters inserted between the sample and the spectrometer block the Rayleigh (elastic) scattered light at the laser wavelength.

The illustration above shows a common system layout in which the Raman emission is collected along a separate optical path from the laser excitation path. Systems designed for imaging (e.g., Raman microscopy systems) or with remote fiber probes are often laid out with the excitation and emission paths coincident, so that both may take advantage of the same fiber and lenses (see Technical Note on page R-9).


There are three basic types of filters used in systems with separate excitation and emission paths: Laser-line filters, Edge Filters, and Notch Filters. The examples, above right, show how the various filters are used. In these graphs the blue lines represent the filter transmission spectra, the green lines represent the laser spectrum, and the red lines

### Laser-line Filter

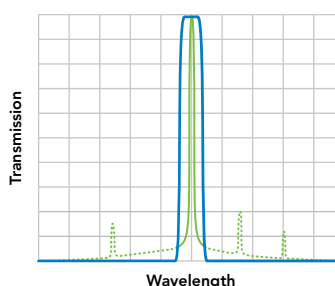


Laser-transmitting filter for both Stokes and Anti-Stokes measurements.

### LWP Edge Filter

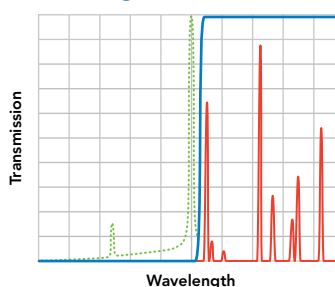

Laser-blocking steep edge filter for superior Stokes measurements.

### Notch Filter

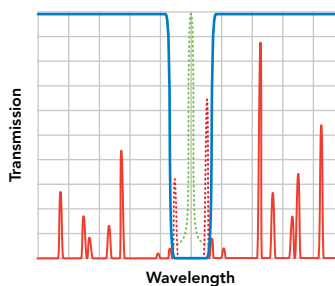

Versatile laser-blocking notch filter for both Stokes and Anti-Stokes measurements.

represent the Raman signal (not to scale).

Laser-line Filters are ideal for use as Laser Transmitting Filters, and Notch Filters are an obvious choice for Laser Blocking Filters. In systems using these two filter types, both Stokes and Anti-Stokes Raman scattering can be measured simultaneously. However, in many cases Edge Filters provide a superior alternative to notch filters. For example, a long-wave-pass (LWP) Edge Filter used as a Laser Blocking Filter for measuring Stokes scattering offers better transmission, higher laser-line blocking, and the steepest edge performance to see Raman signals extremely close to the laser line. For more details on choosing between edge filters and notch filters, see the Technical Note "Edge Filters vs. Notch Filters for Raman Instrumentation" on page R-14.

In systems with a common excitation and emission path, the laser must be introduced into the path with an optic that also allows the Raman emission to be transmitted to the detection system. A 45° dichroic beamsplitter is needed in this case.

If this beamsplitter is not as steep as the edge filter or laser-line filter, one loses the ability to get as close to the laser line as those filters allow.


Semrock manufactures high-performance MaxLine® Laser-line filters (page R-4), RazorEdge® long-wave-pass and short-wave-pass filters (page R-10), EdgeBasic™ value long-wave-pass filters (page R-13), ultrasteep RazorEdge Dichroic™ beamsplitter filters (page R-7), and StopLine® notch filters as standard catalog products. Non-standard wavelengths and specifications for these filters are routinely manufactured for volume OEM applications.

# New – Steeper E-Grade Now Available for MaxLine® Laser-line Filters

## Visible


Wavelength	Guaranteed Transmission	FWHM (nominal)	OD 5 Blue Range (nm)	OD 6 Blue Range (nm)	OD 6 Red Range (nm)	OD 5 Red Range (nm)	12.5 mm Diameter Part Number	25 mm Diameter Part Number
488.0 nm	> 90%	1.4 – 2.0 nm	415.5-449.0	449.0-485.1	490.9-536.8	536.8-625.3	LL01-488E-12.5	LL01-488E-25
532.0 nm	> 90%	1.4 – 2.0 nm	446.5-489.4	489.4-528.8	535.2-585.2	585.2-699.4	LL01-532E-12.5	LL01-532E-25
633.0 nm	> 90%	1.7 – 2.5 nm	515.4-582.2	582.2-629.0	636.6-696.1	696.1-884.7	LL01-633E-12.5	LL01-633E-25
785.0 nm	> 90%	2.3 – 3.5 nm	612.0-722.2	722.2-780.3	789.7-863.5	863.5-1213.8	LL01-785E-12.5	LL01-785E-25

The Figure below illustrates why the MaxLine (LL series, red trace) and RazorEdge (LP series, blue trace) U-grade filters make an ideal filter pair for Raman spectroscopy. The MaxLine filter spectrally “cleans up” the excitation laser light before it reaches the sample under test, allowing only the desired laser line to reach the sample. The RazorEdge filter then removes the laser line from the light scattered by the sample, while efficiently transmitting desired light at wavelengths very close to the laser line. Note the OD 4 separation between the red trace (MaxLine) and the blue trace (U-grade RazorEdge).

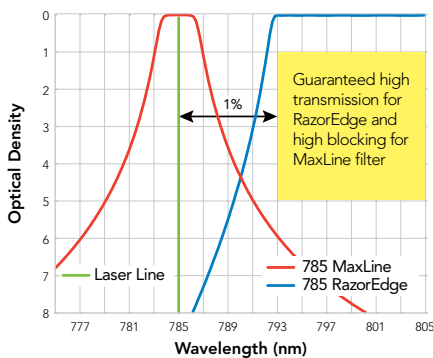

MaxLine and U-grade RazorEdge spectra.


However, this level of cooperativity does not take place between the standard (pre-E) series MaxLine filters and the E-grade RazorEdge filters, which are much steeper than the U-grade filters; see the Figure below, where the LL01-785 cannot adequately clean up the laser light to the red of 785 nm. Note that the LL01-785 has nominal Edge Steepness of 0.9%.

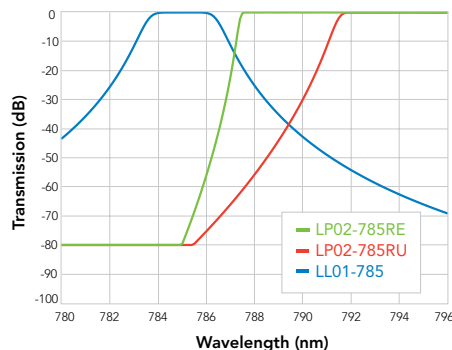

Spectra of MaxLine with U-grade and E-grade RazorEdge filters


The new, steeper MaxLine E-Grade filter corrects this situation due to its increased Edge Steepness (0.2% nominal), allowing better separation from the E-grade RazorEdge.

This new filter subfamily has part numbers with the letter E after the wavelength. In the Figure to the right, the LL01-785E is distinguished from the standard steepness LL01-785 shown in the Figure below.

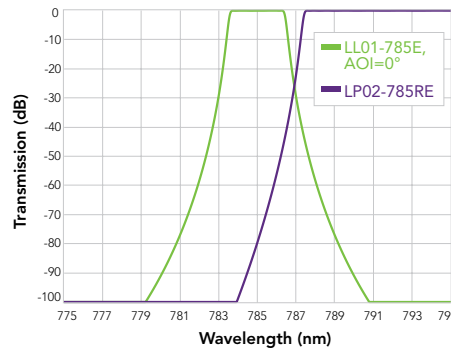

Superposed graphs of E-grade MaxLine and RazorEdge spectra.

The built-in angle-tilt feature is the other novel feature of the MaxLine E-grade series. With reference to the Figure below, the separation between the RazorEdge and MaxLine filters can be increased by adjusting the AOI of the MaxLine filter to shift the steep red edge to a lower wavelength. This is shown in the Figure below, where about 4 OD of suppression is achieved at the crossover wavelength at the typical maximum rotation of 5°.

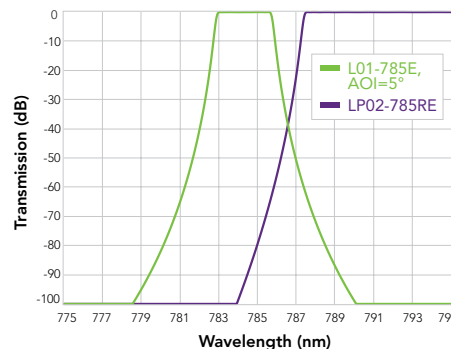


Shift of theoretical spectra of LL01-785E with AOI

# MaxLine® Laser-line Filters

Semrock MaxLine Laser-line Filters have an unprecedented high transmission exceeding 90% at the laser line, while rapidly rolling off to an optical density (OD) > 5 at wavelengths differing by only 1% from the laser wavelength, and OD > 6 at wavelengths differing by only 1.5% from the laser wavelength. U.S. Patent No. 7,119,960.

- Highest laser-line transmission – stop wasting expensive laser light
- Steepest edges – perfect match to RazorEdge® U-grade filters (see page R-13)
- Ideal complement to StopLine® deep notch filters for fluorescence and other applications
- Hard dielectric coatings for proven reliability and durability
- For diode lasers, use our MaxDiode™ Laser Clean-up filters

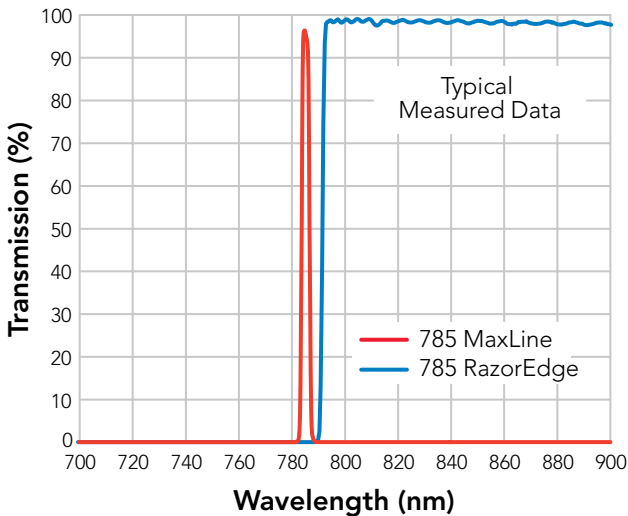

## ★ PRODUCT NOTE


### RazorEdge and MaxLine® are a Perfect Match

The MaxLine (see page R-3) and RazorEdge U-grade (see page R-14) filters make an ideal filter pair for applications like Raman spectroscopy – they fit together like hand-in-glove. The MaxLine filter spectrally “cleans up” the excitation laser light before it reaches the sample under test – allowing only the desired laser line to reach the sample – and then the RazorEdge filter removes the laser line from the light scattered off the sample, while efficiently transmitting desired light at wavelengths very close to the laser line.

Typical measured spectral curves of 785 nm filters on a linear transmission plot demonstrate how the incredibly steep edges and high transmission exhibited by both of these filters allow them to be spectrally positioned very close together, while still maintaining complementary transmission and blocking characteristics.

The optical density plot illustrates the complementary nature of these filters on a logarithmic scale using the theoretical design spectral curves. The MaxLine filter provides very high transmission (> 90%) of light immediately in the vicinity of the laser line, and then rapidly rolls off to achieve very high blocking (> OD 5) at wavelengths within 1% of the laser line. The RazorEdge filter provides extremely high blocking (> OD 6) of the laser line itself, and then rapidly climbs to achieve very high transmission (> 90%) of the desired signal light at wavelengths only 1% away from the laser line.

If you are currently using an E-grade RazorEdge filter and need a laser clean-up filter, see our new E-grade Maxline filters, page R-3.

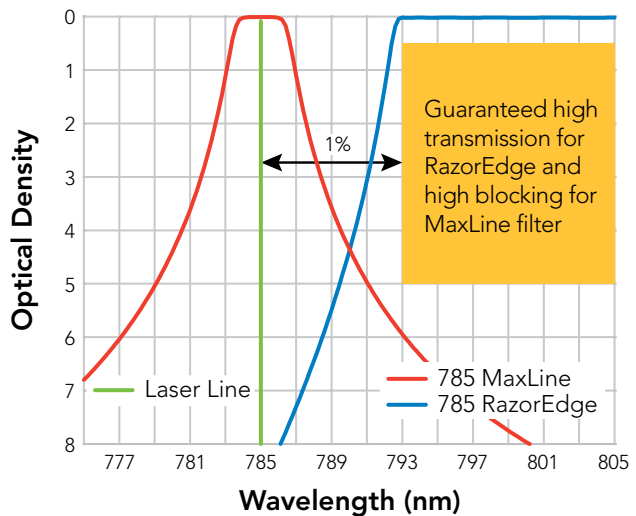

# MaxLine® Laser-line Filters (Continued)

## Ultraviolet

Wavelength	Guaranteed Transmission	Typical Bandwidth	OD 5 Blue Range (nm)	OD 6 Blue Range (nm)	OD 6 Red Range (nm)	OD 5 Red Range (nm)	12.5 mm Diameter Part Number	25 mm Diameter Part Number
248.6 nm	> 40%	1.7 nm	228.2-246.1	228.7-244.9	252.3-273.5	251.1-279.9	LL01-248-12.5	LL01-248-25
266.0 nm	> 55%	1.9 nm	242.8-263.3	244.7-262.0	270.0-292.6	268.7-302.2	LL01-266-12.5	LL01-266-25
280.0 nm	> 60%	1.5 nm	254.4-277.2	257.6-275.8	282.8-320.4	284.2-308	LL01-280-12.5	LL01-280-25
320.0 nm	> 70%	1.3 nm	286.9-316.8	294.4-315.2	323.2-373.8	324.8-352	LL01-320-12.5	LL01-320-25
325.0 nm	> 80%	1.2 nm	291.0-321.8	299.0-320.1	329.9-357.5	328.3-380.7	LL01-325-12.5 ★	LL01-325-25
355.0 nm	> 80%	1.3 nm	314.8-351.5	326.6-349.7	360.3-390.5	358.6-422.5	LL01-355-12.5	LL01-355-25
360.0 nm	> 85%	1.1 nm	318.7-356.4	331.2-354.6	363.6-429.6	365.4-396	LL01-360-12.5	LL01-360-25

## Visible

Wavelength	Guaranteed Transmission	Typical Bandwidth	OD 5 Blue Range (nm)	OD 6 Blue Range (nm)	OD 6 Red Range (nm)	OD 5 Red Range (nm)	12.5 mm Diameter Part Number	25 mm Diameter Part Number
405.0 nm	> 90%	1.5 nm	353.5-401.0	372.6-398.9	411.1-445.5	409.1-495.3	LL01-405-12.5 ★	LL01-405-25
407.0 nm	> 90%	1.5 nm	355.0-402.9	374.4-400.9	413.1-447.7	411.1-498.3	LL01-407-12.5	LL01-407-25
441.6 nm	> 90%	1.7 nm	381.0-437.2	406.3-435.0	448.2-485.8	446.0-551.1	LL01-442-12.5	LL01-442-25
457.9 nm	> 90%	1.7 nm	393.1-453.3	421.3-451.0	464.8-503.7	462.5-576.7	LL01-458-12.5	LL01-458-25
473.0 nm	> 90%	1.8 nm	404.2-468.3	435.2-465.9	480.1-520.3	477.7-600.9	LL01-473-12.5	LL01-473-25
488.0 nm	> 90%	1.9 nm	415.1-483.1	449.0-480.7	495.3-536.8	492.9-625.3	LL01-488-12.5 ★	LL01-488-25
488.0 nm	> 90%	1.4 – 2.0 nm	415.5-449.0	449.0-485.1	490.9-536.8	536.8-625.3	LL01-488E-12.5	LL01-488E-25
514.5 nm	> 90%	2.0 nm	434.1-509.4	473.3-506.8	522.2-566.0	519.6-669.5	LL01-514-12.5	LL01-514-25
532.0 nm	> 90%	2.0 nm	446.5-526.7	489.4-524.0	540.0-585.2	537.3-699.4	LL01-532-12.5 ★	LL01-532-25
532.0 nm	> 90%	1.4 – 2.0 nm	446.5-489.4	489.4-528.8	535.2-585.2	585.2-699.4	LL01-532E-12.5	LL01-532E-25
543.5 nm	> 90%	2.1 nm	454.6-538.1	500.0-535.3	551.7-597.9	548.9-719.5	LL01-543-12.5	LL01-543-25
561.4 nm	> 90%	2.1 nm	467.0-555.8	516.5-553.0	569.8-617.5	567.0-751.2	LL02-561-12.5 ★	LL02-561-25
568.2 nm	> 90%	2.2 nm	471.7-562.5	522.7-559.7	576.7-625.0	573.9-763.4	LL01-568-12.5	LL01-568-25
632.8 nm	> 90%	2.4 nm	515.4-626.5	582.2-623.3	642.3-696.1	639.1-884.7	LL01-633-12.5 ★	LL01-633-25
633.0 nm	> 90%	1.7 – 2.5 nm	515.4-582.2	582.2-629.0	636.6-696.1	696.1-884.7	LL01-633E-12.5	LL01-633E-25
638.0 nm	> 90%	2.4 nm	518.8-631.6	587-628.4	647.6-701.8	644.4-894.9	LL01-638-12.5	LL01-638-25
647.1 nm	> 90%	2.5 nm	524.8-640.6	595.3-637.4	656.8-711.8	653.6-912.9	LL01-647-12.5	LL01-647-25
671.0 nm	> 90%	2.6 nm	540.4-664.3	617.3-660.9	681.1-738.1	677.7-961.2	LL01-671-12.5	LL01-671-25
785.0 nm	> 90%	2.3 – 3.5 nm	612.0-722.2	722.2-780.3	789.7-863.5	863.5-1213.8	LL01-785E-12.5	LL01-785E-25

## Near-Infrared

Wavelength	Guaranteed Transmission	Typical Bandwidth	OD 5 Blue Range (nm)	OD 6 Blue Range (nm)	OD 6 Red Range (nm)	OD 5 Red Range (nm)	12.5 mm Diameter Part Number	25 mm Diameter Part Number
780.0 nm	> 90%	3.0 nm	609.0-772.2	717.6-768.3	791.7-858.0	787.8-1201.8	LL01-780-12.5	LL01-780-25
785.0 nm	> 90%	3.0 nm	612.0-777.2	722.2-773.2	796.8-863.5	792.9-1213.8	LL01-785-12.5 ★	LL01-785-25
808.0 nm	> 90%	3.1 nm	625.9-799.9	743.4-795.9	820.1-888.8	816.1-1033.5	LL01-808-12.5 ★	LL01-808-25
810.0 nm	> 90%	3.1 nm	627.1-801.9	745.2-797.9	822.2-891.0	818.1-1143.4	LL01-810-12.5 ★	LL01-810-25
830.0 nm	> 90%	3.2 nm	639.1-821.7	763.6-817.6	842.5-913.0	838.3-1067.9	LL01-830-12.5	LL01-830-25
852.0 nm	> 90%	3.2 nm	652-843.5	783.8-839.2	864.8-937.2	860.5-1106.6	LL01-852-12.5	LL01-852-25
976.0 nm	> 90%	3.7 nm	722.2-966.2	897.9-961.4	990.6-1073.6	985.8-1325.2	LL01-976-12.5	LL01-976-25
980.0 nm	> 90%	3.7 nm	724.4-970.2	901.6-965.3	994.7-1078.0	989.8-1332.6	LL01-980-12.5	LL01-980-25
1030.0 nm	> 90%	3.9 nm	1014.6-1019.7	947.6-1014.6	1045.5-1133	1040.3-1368.2	LL01-1030-12.5	LL01-1030-25
1047.1 nm	> 90%	4.0 nm	963.3-1036.6	963.3-1031.4	1062.8-1151.8	1057.6-1398.6	LL01-1047-12.5	LL01-1047-25
1064.0 nm	> 90%	4.0 nm	978.9-1053.4	978.9-1048.0	1080.0-1170.4	1074.6-1428.9	LL01-1064-12.5 ★	LL01-1064-25

Price Visit [index-hs.com/semrock](http://index-hs.com/semrock) to Search the MaxLine Filter Family

★ = High Selling Item – quantity ≤ 3 usually ships in 3 days, demand dependent.

## ★ PRODUCT NOTE

# Edge Steepness and Transition Width

Semrock edge filters, including our steepest RazorEdge® Raman filters as well as our EdgeBasic™ filters for application-specific Raman systems and fluorescence imaging, are specified with a guaranteed “Transition Width.”

**Transition Width** = maximum allowed spectral width between the laser line (where OD > 6) and the 50% transmission point

Any given filter can also be described by its “Edge Steepness,” which is the actual steepness of the filter, regardless of the precise wavelength placement of the edge.

**Edge Steepness** = actual steepness of a filter measured from the OD 6 point to the 50% transmission point

Figure 1 illustrates Transition Width and Edge Steepness for an edge filter designed to block the 785 nm laser line (example shows a U-grade RazorEdge filter). The table below lists the guaranteed Transition Width and typical Edge Steepness (for 25 mm diameter parts) for Semrock edge filters.

All Verona and RazorEdge filters provide exceptional steepness to allow measurement of signals very close to the blocked laser line with high signal-to-noise ratio. However, the state-of-the-art E-grade RazorEdge filters take closeness to an extreme level.

Figure 2 illustrates that U-grade RazorEdge filters have a transition width that is 1% of the laser wavelength. E-grade filters have a transition width that is twice as narrow, or 0.5% of the laser line!

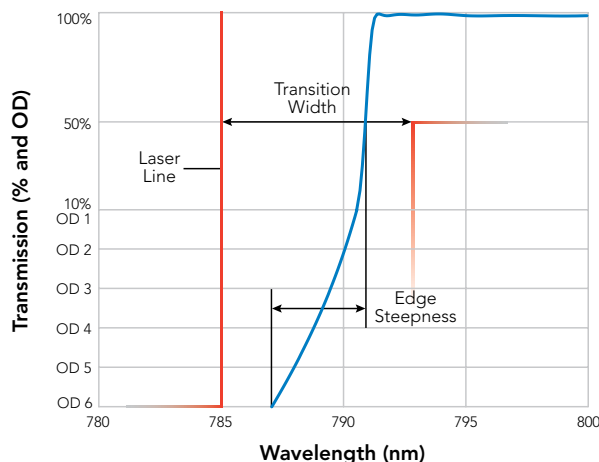


Figure 1: Transition width and edge steepness illustrated.

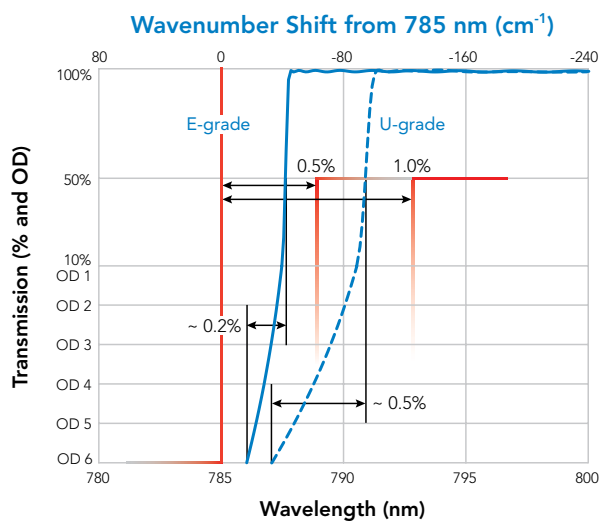


Figure 2: Transition widths and edge steepnesses for LP02-785RE and LP02-785RU filters (see page R-11).

Edge Filter Type	Guaranteed Transition Width (% of laser wavelength)	Typical Edge Steepness (% of laser wavelength)
Verona	Refer to page R-8	
RazorEdge E-grade	< 0.5% (< 90 cm <sup>-1</sup> for 532)	0.2% (1.1 nm for 532)
RazorEdge U-grade	< 1.0% (< 186 cm <sup>-1</sup> for 532)	0.5% (2.7 nm for 532)
EdgeBasic	< 2.5% (< 458 cm <sup>-1</sup> for 532)	1.5% (8.0 nm for 532)
* except UV filters		



# RazorEdge® Dichroic™ Beamsplitters

The unique RazorEdge Dichroic beamsplitters exhibit unparalleled performance. Each filter reflects a standard laser line incident at 45° while efficiently passing the longer Raman-shifted wavelengths. They exhibit ultrasteep transition from reflection to transmission, far superior to anything else available on the open market. The guaranteed transition width of < 1% of the laser wavelength for U-grade (regardless of polarization) makes these filters a perfect match to our popular normal-incidence RazorEdge ultrasteep long-wave-pass filters. These beamsplitters are so innovative that they are patent pending.



Available as Either Mounted in 25 mm Diameter x 3.5 mm Thick Black-Anodized Aluminum Ring or Unmounted as 25.2 x 35.6 x 1.1 mm or 25.2 x 35.6 x 2.0 mm

Laser Line	Transition Width	Passband	25 mm Mounted Part Number	25.2 x 35.6 x 1.1 mm Part Number	25.2 x 35.6 x 2.0 mm Part Number
532.0 nm	< 186 cm <sup>-1</sup>	538.9 – 824.8 nm	LPD02-532RU-25	LPD02-532RU-25x36x1.1 ★	LPD02-532RU-25x36x2.0
632.8 nm	< 156 cm <sup>-1</sup>	641.0 – 980.8 nm	LPD02-633RU-25	LPD02-633RU-25x36x1.1 ★	LPD02-633RU-25x36x2.0
785.0 nm	< 126 cm <sup>-1</sup>	795.2 – 1213.8 nm	LPD02-785RU-25	LPD02-785RU-25x36x1.1 ★	LPD02-785RU-25x36x2.0
830.0 nm	< 119 cm <sup>-1</sup>	840.8 – 1286.5 nm	LPD02-830RU-25	LPD02-830RU-25x36x1.1	LPD02-830RU-25x36x2.0
1064.0 nm	< 93 cm <sup>-1</sup>	1077.8 – 1650.8 nm	LPD02-1064RU-25		LPD02-1064RU-25x36x2.0

★ = High Selling Item – quantity ≤ 3 usually ships in 3 days, demand dependent.

 AVAILABLE IN 1.1 mm THICKNESSES FOR MICROSCOPES

 SEE SPECTRAL GRAPHS AND GET NUMERIC DATA FOR ALL OUR FILTERS AT [IDEX-HS.COM/SEMROCK](http://IDEX-HS.COM/SEMROCK)

## Dichroic Beamsplitter Specifications



Property	Specification	Comment	
Edge Steepness (typical)	0.5% of laser wavelength (2.5 nm or 90 cm <sup>-1</sup> for 532 nm filter)	Measured from 5% to 50% transmission for light with average polarization	
Transition Width	U-grade < 1% of laser wavelength	Measured from laser wavelength to 50% transmission wavelength for light with average polarization	
Reflection at Laser	Wavelength > 98% (s-polarization) > 90% (p-polarization)		
Average Passband	Transmission > 93%	Averaged over the Passband (Passband wavelengths detailed above)	
Dependence of Wavelength on Angle of Incidence (Edge Shift)	≤ 0.2% / degree	Linear relationship valid between 35° & 55° (see MyLight for actual performance)	
Cone Half Angle (for non-collimated light)	< 0.5°	Rays uniformly distributed and centered at 45°	
Size of Round Dichroics	Clear Aperture	≥ 22 mm	
	Outer Diameter	25.0 + 0.0 / - 0.1 mm	Black-anodized aluminum ring
	Overall Thickness	3.5 ± 0.1 mm	Black-anodized aluminum ring
Size of Rectangular Dichroics	Clear Aperture	> 80%	Elliptical
	Size	25.2 mm x 35.6 mm ± 0.1 mm	
	Unmounted Thickness	1.05 mm ± 0.05 mm	
Wedge Angle	< 20 arcseconds		
Flatness	Reflection of a collimated, Gaussian laser beam with waist diameter up to 3 mm causes less than one Rayleigh Range of focal shift after a focusing lens.		

# Verona™ Raman Filter Family

Semrock's long-pass Verona optical filters provide industry-best solutions for deep blocking close to the laser line for Raman applications. Verona is just the beginning of our evolving Raman product line to match the rapidly developing Raman market. Our understanding of how crucial flat transmission is for the maximal collection of weak Raman signals is why we offer a filter with both low ripple and improved steepness from our RazorEdge™ series, making the Verona family the new gold standard for Raman analysis.

Push the limits of what you can see with your Raman system with Verona Raman filters.

- Transition width  $\leq 0.2\%$  relative to the edge wavelength while still achieving transmission  $> 90\%$  and blocking  $> OD6$
- Steepness significantly improved over our RazorEdge series
- Low ripple to provide the best signal-to-noise ratio and allow for maximal collection of weak Raman spectral features
- 785 nm, 633 nm, 532 nm, and 488 nm wavelengths are now available
- Deep blocking of  $> OD6$  at the laser line that eliminates bleed-through of excitation light
- Hard coated for high laser damage threshold to reduce degradation in performance under ambient conditions
- With our award-winning, proprietary KolaDeep™ spectral measurement system, we provide the spectral data to prove our superior deep blocking and edge steepness performance
- Available in 12.5 mm housed standard sized parts to meet your application needs

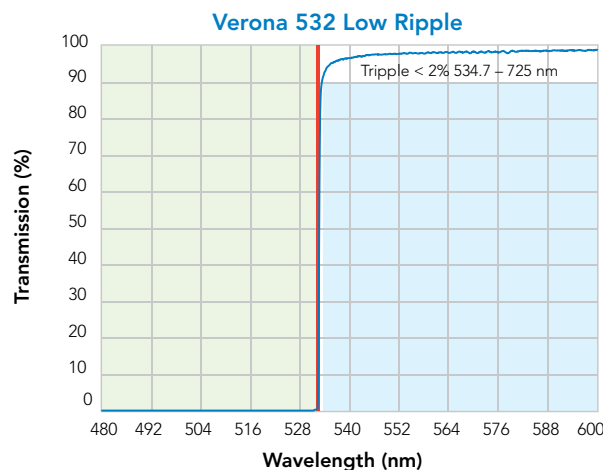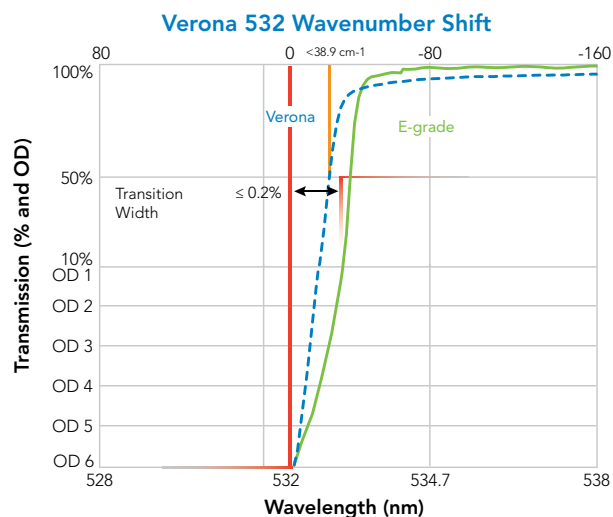

## 12.5 mm Diameter

Laser Line	Transition Width	Passband	Part Number	Price
488 nm	$< 38.8 \text{ cm}^{-1}$	490.4 – 1100.8 nm	VLP01-488-12.5	\$815
532 nm	$< 38.9 \text{ cm}^{-1}$	534.7 – 1300 nm	VLP01-532-12.5	\$815
633 nm	$< 25.9 \text{ cm}^{-1}$	636.2 – 1427.4 nm	VLP02-633-12.5	\$815
785 nm	$< 25.9 \text{ cm}^{-1}$	788.9 – 1770.7 nm	VLP02-785-12.5	\$815

## Verona Specifications

Property	Specifications	Comments
Edge Steepness (Typical)	0.2% of laser wavelength	Measured from OD = 6 to 50% transmission wavelength
Ripple on passband Transmission	$< 2\%$	
Angle of Incidence	$0.0^\circ \pm 2.0^\circ$	Range for above optical specifications
Cone Half Angle	$0^\circ$	
Clear Aperture	$\geq 10 \text{ mm}$	
Outer Diameter	$12.5 \text{ mm} + 0.0 / - 0.1 \text{ mm}$	
Substrate Thickness	3.0 mm	
Mounted Thickness	5.0 mm	



# RazorEdge® Common Specifications

## RazorEdge Specifications

Properties apply to all long-wave-pass and short-wave-pass edge filters unless otherwise noted

Property	Specification	Comment
Edge Steepness (typical)	E-grade 0.2% of laser wavelength U-grade 0.5% of laser wavelength	Measured from OD 6 to 50%; Up to 0.8% for 248-300 nm filters and 3.3% for 224 nm filter
Blocking at Laser Wavelength	> 6 OD	OD = -log <sub>10</sub> (transmission)
Transition Width	E-grade < 0.5% of laser wavelength U-grade < 1% of laser wavelength	Measured from laser wavelength to 50% transmission wavelength; < 4.5% for 224 nm filter
Guaranteed Passband Transmission	> 93%	Except > 90% for 224 – 325 nm filters; Averaged over the Passband
Typical Passband Transmission	> 98%	
Angle of Incidence	0.0° ± 2.0°	Range for above optical specifications
Cone Half Angle	< 5°	Rays uniformly distributed about 0°
Angle Tuning Range <sup>1</sup>	-0.3% of Laser Wavelength	Wavelength "blue shift" attained by increasing angle from 0° to 8°
Laser Damage Threshold	0.5 J/cm <sup>2</sup> @ 266 nm 1 J/cm <sup>2</sup> @ 532 nm	10 ns pulse width Tested for 266 and 532 nm filters only (see page 109)
Clear Aperture	> 22 mm (or > 45 mm)	
Outer Diameter	25.0 + 0.0 / -0.1 mm (or 50.0 + 0.0 / -0.1 mm)	Black-anodized aluminum ring
Substrate Thickness	2.0 mm	
Overall Thickness	3.5 ± 0.1 mm	Black-anodized aluminum ring (thickness measured unmounted)
Beam Deviation	< 10 arcseconds	

<sup>1</sup> For small angles (in degrees), the wavelength shift near the laser wavelength is  $\Delta\lambda$  (nm) =  $-5.0 \times 10^{-5} \times \lambda \times q^2$  and the wavenumber shift is  $\Delta(\text{wavenumbers})$  (cm<sup>-1</sup>) =  $500 \times q^2 / \lambda$ , where  $\lambda$  (in nm) is the laser wavelength. See Wavenumbers Technical Note on page 102.

## General Specifications (all RazorEdge filters)


Property	Specification	Comment
Coating Type	"Hard" ion-beam-sputtered	
Reliability and Durability	Ion-beam-sputtered, hard-coated technology with epoxy-free, single-substrate construction for unrivaled filter life. RazorEdge filters are rigorously tested and proven to MIL-STD-810F and MIL-C-48497A environmental standards.	
Transmitted Wavefront Error	< $\lambda / 4$ RMS at $\lambda = 633$ nm	Peak-to-valley error < 5 x RMS value measured within clear aperture
Surface Quality	60-40 scratch-dig	
Temperature Dependence	< 5 ppm / °C	
Substrate Material	Ultra-low autofluorescence fused silica	
Filter Orientation	For mounted filters, arrow on ring indicates preferred direction of propagation of transmitted light. For rectangular dichroics, reflective coating side should face toward light source and sample.	

## TECHNICAL NOTE

# RazorEdge Filter Layouts

Only the unique RazorEdge Dichroic beamsplitter reflects a standard laser line incident at 45° while transmitting longer Raman-shifted wavelengths with an ultrasteep transition far superior to anything else available on the open market. The guaranteed transition width of < 1% of the laser wavelength for U-grade (regardless of polarization) makes these filters a perfect match to our popular normal-incidence RazorEdge ultrasteep long-wave-pass filters.

In order for the two-filter configuration to work, the 45° beamsplitter must be as steep as the laser-blocking filter. Traditionally thin-film filters could not achieve very steep edges at 45° because of the "polarization splitting" problem – the edge position tends to be different for different polarizations of light. However, through continued innovation in thin-film filter technology, Semrock has been able to achieve ultrasteep 45° beamsplitters with the same steepness of our renowned RazorEdge laser-blocking filters: the transition from the laser line to the passband of the filter is guaranteed to be less than 1% of the laser wavelength (for U-grade filters).


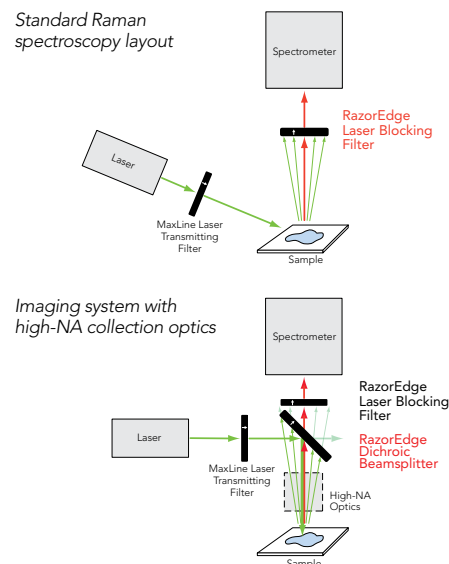


# RazorEdge® Long Wave Pass Raman Edge Filters

Semrock stocks an unsurpassed selection of the highest performance edge filters available for Raman Spectroscopy, with edge wavelengths from 224 to 1319 nm. Now you can see the weakest signals closer to the laser line than you ever have before. With their deep laser-line blocking, ultra-wide and low-ripple passbands, proven hard-coating reliability, and high laser damage threshold, they offer performance that lasts. U.S. Patent No. 7,068,430.

- The steepest edge filters on the market – RazorEdge E-grade filters, see how steep on page R-4
- For long-wave-pass edge filters and normal incidence, see right
- For short-wave-pass edge filters and normal incidence, see page R-11
- For ultrasteep 45° beamsplitters, see page R-7
- For a suitably matched laser-line filter, see page R-3 and R-4

## ACTUAL MEASURED DATA RAZOREDGE RAMAN FILTER SPECTRA

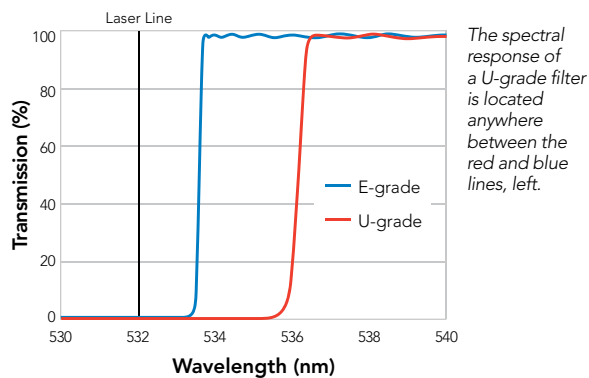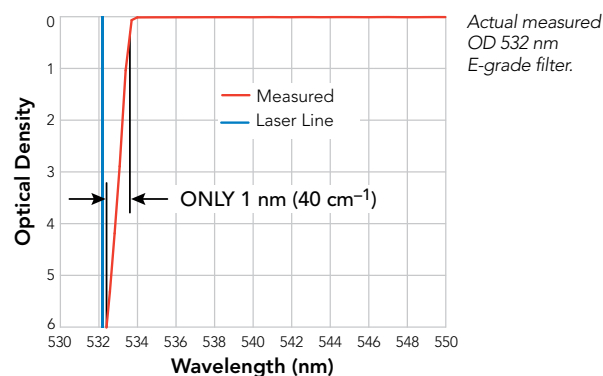

➤ FOR 50 MM LWP EDGE FILTERS, REFER TO EACH FILTER'S PRODUCT PAGE AT [WWW.IDEX-HS.COM/SEMROCK](http://WWW.IDEX-HS.COM/SEMROCK)

Note: E-Grade LWP filters are available in sizes 10 mm to 25 mm.

### 25 mm Diameter

Laser Line	Transition Width <sup>1</sup>	Passband	Part Number	Price
224.3 nm	< 1920 cm <sup>-1</sup>	235.0-505.9 nm	LP02-224R-25	\$1,230
244 nm	< 498 cm <sup>-1</sup>	247.6-550.4 nm	LP02-244RS-25	\$1,230
248.6 nm	< 805 cm <sup>-1</sup>	261.0-560.8 nm	LP02-248RS-25	\$1,230
266.0 nm	< 372 cm <sup>-1</sup>	272.4-600.0 nm	LP02-266RU-25	\$1,095
280.0 nm	< 354 cm <sup>-1</sup>	283.6 – 631.6 nm	LP02-280RU-25	\$745
320.0 nm	< 309 cm <sup>-1</sup>	324.2 – 721.8 nm	LP02-320RU-25	\$745
325.0 nm	< 153 cm <sup>-1</sup> < 305 cm <sup>-1</sup>	327.1-733.1 nm 329.2-733.1 nm	LP03-325RE-25 LP03-325RU-25	\$1,095 \$745
355.0 nm	< 140 cm <sup>-1</sup> < 279 cm <sup>-1</sup>	357.3-800.8 nm 359.6-800.8 nm	LP02-355RE-25 LP02-355RU-25	\$1,230 \$745
360.0 nm	< 275 cm <sup>-1</sup>	364.7 – 812 nm	LP02-360RU-25	\$745
405.0 nm	< 244 cm <sup>-1</sup>	410.3 – 913.5 nm	LP02-405RU-25	\$745
407.0 nm	< 243 cm <sup>-1</sup>	412.3-918.0 nm	LP02-407RU-25	\$745
441.6 nm	< 224 cm <sup>-1</sup>	447.3-996.1 nm	LP02-442RU-25	\$850
457.9 nm	< 109 cm <sup>-1</sup> < 216 cm <sup>-1</sup>	460.9-1032.9 nm 463.9-1032.9 nm	LP03-458RE-25 LP03-458RU-25	\$1,230 \$745
473.0 nm	< 105 cm <sup>-1</sup> < 209 cm <sup>-1</sup>	476.1-1066.9 nm 479.1-1066.9 nm	LP02-473RE-25 LP02-473RU-25	\$1,095 \$745
488.0 nm	< 102 cm <sup>-1</sup> < 203 cm <sup>-1</sup>	491.2-1100.8 nm 494.3-1100.8 nm	LP02-488RE-25 LP02-488RU-25	\$1,095 \$745 ★
514.5 nm	< 97 cm <sup>-1</sup> < 192 cm <sup>-1</sup>	517.8-1160.5 nm 521.2-1160.5 nm	LP02-514RE-25 LP02-514RU-25	\$1,230 \$850
532.0 nm	< 90 cm <sup>-1</sup> < 186 cm <sup>-1</sup>	535.4-1200.0 nm 538.9-1200.0 nm	LP03-532RE-25 LP03-532RU-25	\$1,095 \$745 ★
561.4 nm	< 89 cm <sup>-1</sup> < 176 cm <sup>-1</sup>	565.0-1266.3 nm 568.7-1266.3 nm	LP02-561RE-25 LP02-561RU-25	\$1,230 \$850
632.8 nm	< 79 cm <sup>-1</sup> < 156 cm <sup>-1</sup>	636.9-1427.4 nm 641.0-1427.4 nm	LP02-633RE-25 LP02-633RU-25	\$1,095 \$745 ★
638 nm	< 78 cm <sup>-1</sup> < 155 cm <sup>-1</sup>	642.1-1439.1 nm 646.3-1439.1 nm	LP02-638RE-25 LP02-638RU-25	\$1,230 \$850
647.1 nm	< 153 cm <sup>-1</sup>	655.5-1459.6 nm	LP02-647RU-25	\$850
664.0 nm	< 149 cm <sup>-1</sup>	672.6-1497.7 nm	LP02-664RU-25	\$850
671.0 nm	< 147.6 cm <sup>-1</sup>	679.7-1513.5 nm	LP02-671RU-25	\$745
780.0 nm	< 127 cm <sup>-1</sup>	790.1-1759.4 nm	LP02-780RU-25	\$850
785.0 nm	< 63 cm <sup>-1</sup> < 126 cm <sup>-1</sup>	790.1-1770.7 nm 795.2-1770.7 nm	LP02-785RE-25 LP02-785RU-25	\$1,095 \$745 ★
808.0 nm	< 62 cm <sup>-1</sup> < 123 cm <sup>-1</sup>	813.3-1822.6 nm 818.5-1822.6 nm	LP02-808RE-25 LP02-808RU-25	\$1,230 \$850
830.0 nm	< 60 cm <sup>-1</sup> < 119 cm <sup>-1</sup>	835.4-1872.2 nm 840.8-1872.2 nm	LP02-830RE-25 LP02-830RU-25	\$1,230 \$745
980.0 nm	< 51 cm <sup>-1</sup> < 101 cm <sup>-1</sup>	986.4-2000.0 nm 992.7-2000.0 nm	LP02-980RE-25 LP02-980RU-25	\$1,230 \$745
1064.0 nm	< 47 cm <sup>-1</sup> < 93 cm <sup>-1</sup>	1070.9-2000.0 nm 1077.8-2000.0 nm	LP02-1064RE-25 LP02-1064RU-25	\$1,230 \$745 ★


★ = High Selling Item – quantity ≤ 3 usually ships in 3 days, demand dependent.

# RazorEdge® Short Wave Pass Raman Edge Filters

These unique filters are ideal for Anti-Stokes Raman applications. An addition to the popular high-performance RazorEdge family of steep edge filters, these short-wave-pass filters are designed to attenuate a designated laser-line by six orders of magnitude, and yet maintain a typical edge steepness of only 0.5% of the laser wavelength. Both short and long-wave-pass RazorEdge filters are perfectly matched to Semrock's popular MaxLine® laser-line cleanup filters. U.S. Patent No. 7,068,430



## ACTUAL MEASURED DATA FROM A 632.8 nm RAZOREEDGE FILTER

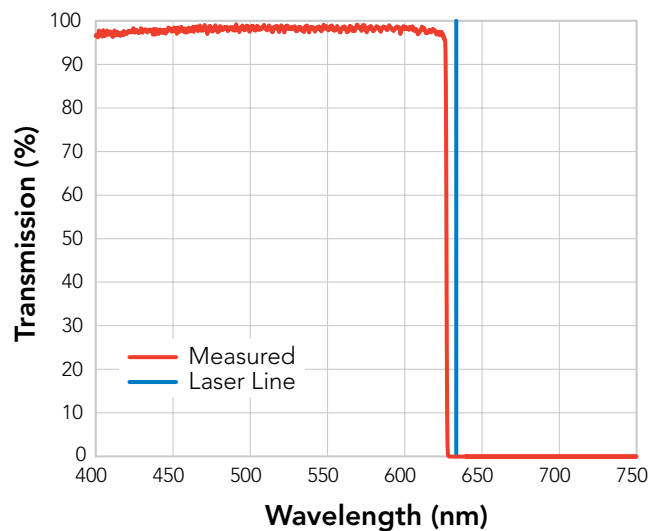


SEE SPECTRAL GRAPHS AND GET NUMERIC DATA FOR ALL OUR FILTERS AT [IDEX-HS.COM/SEMROCK](http://IDEX-HS.COM/SEMROCK)

### 25 mm Diameter

Laser Line	Transition Width	Passband	Part Number	Price
532.0 nm	< 186 cm <sup>-1</sup>	350.0 – 525.2 nm	SP01-532RU-25	\$855
561.4 nm	< 176 cm <sup>-1</sup>	400.0 – 554.1 nm	SP01-561RU-25	\$745
632.8 nm	< 160 cm <sup>-1</sup>	372.0 – 624.6 nm	SP01-633RU-25	\$855
785.0 nm	< 129 cm <sup>-1</sup>	400.0 – 774.8 nm	SP01-785RU-25	\$745

FOR 50 MM LWP EDGE FILTERS, REFER TO EACH FILTER'S PRODUCT PAGE AT [WWW.IDEX-HS.COM/SEMROCK](http://WWW.IDEX-HS.COM/SEMROCK)

# EdgeBasic™ Long/Short Wave Pass Filters Common Specifications

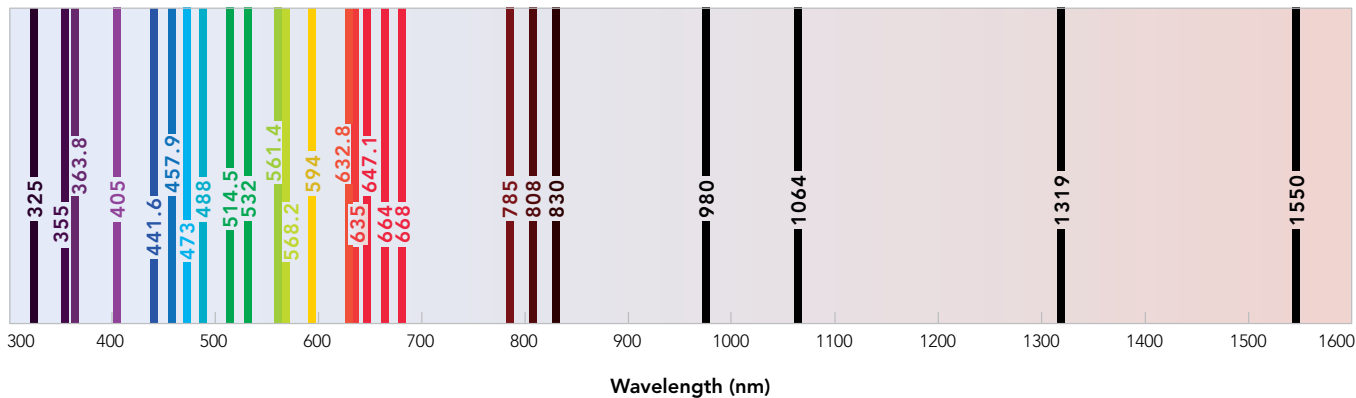

## Longpass specifications

Property	Value	Comments
Edge Steepness (typical)	1.5% of $\lambda_{\text{long}}$	Measured from OD 6 to 50%
Transition Width	< 2.5% of $\lambda_{\text{long}}$	From $\lambda_{\text{long}}$ to the 50% transmission wavelength
Blocking at Laser Wavelengths	ODabs > 6 from 80% of $\lambda_{\text{short}}$ to $\lambda_{\text{long}}$ ODavg > 5 from 270 nm to 80% of $\lambda_{\text{short}}$ ( $\lambda_s \leq 1064$ nm) ODavg > 5 from 800 nm to 80% of $\lambda_{\text{short}}$ ( $\lambda_s > 1064$ nm)	OD = - log <sub>10</sub> (transmission)
Guaranteed Transmission	> 93%	Averaged over the passband

## Shortpass Specifications

Property	Value	Comments
Edge Steepness (typical)	1.5% of $\lambda_{\text{short}}$	Measured from OD 6 to 50%
Transition Width	< 2.5% of $\lambda_{\text{short}}$	From 50% transmission wavelength to $\lambda_{\text{short}}$
Blocking at Laser Wavelengths	ODabs > 6 from $\lambda_{\text{short}}$ to 120% of $\lambda_{\text{long}}$ ODavg > 5 from 120% of $\lambda_{\text{long}}$ to 750 nm ODavg > 4 from 750 nm to 925 nm ODavg > 3 from 925 nm to 1200 nm	OD = - log <sub>10</sub> (transmission)
Guaranteed Transmission	> 93%	Averaged over the passband >400 nm

## Common Specifications

Property	Value	Comments
Guaranteed Transmission	> 93%	Averaged over the passband for Shortpass > 80% 350 – 400 nm
Typical Transmission	> 98%	Averaged over the passband
Angle of Incidence	0.0° ± 2.0°	Range for above optical specifications
Cone Half Angle	< 5°	Rays uniformly distributed about 0°
Angle Tuning Range	~ 0.3% of Laser Wavelength	Wavelength “blue shift” increasing angle from 0° to 8°
Substrate Material	Low-autofluorescence optical quality glass	
Substrate Thickness	2.0 ± 0.1 mm	
Clear Aperture	> 22 mm	
Outer Diameter	25.0 + 0.0 / - 0.1 mm	Black-anodized aluminum ring
Overall Thickness	3.5 ± 0.1 mm	Black-anodized aluminum ring
Beam Deviation	< 10 arc seconds	
Surface Quality	60-40 scratch-dig	
Filter Orientation	Arrow on ring indicates preferred direction of propagation of light	

# EdgeBasic™ Long / Short Wave Pass Filters

EdgeBasic long wave pass and short wave pass filters offer a superb combination of performance and value for applications in Raman spectroscopy and fluorescence imaging and measurements. This group of filters is ideal for specific Raman applications that do not require measuring the smallest possible Raman shifts, yet demand exceptional laser-line blocking and high transmission over a range of Raman lines.

- › Deep laser-line blocking – for maximum laser rejection (OD > 6)
- › Extended short-wavelength blocking (LWP) – for high-fidelity fluorescence imaging
- › High signal transmission – to detect the weakest signals (> 98% typical)
- › Proven no burn-out durability – for lasting and reliable performance
- › For the ultimate performance, upgrade to state-of-the-art RazorEdge® Raman filters



## EdgeBasic Long Wave Pass

Nominal Laser Wavelength	Laser Wavelength Range		Passband	Part Number	Price
	$\lambda$ short	$\lambda$ long			
325 nm	325.0 nm	325.0 nm	334.1 – 900.0 nm	BLP01-325R-25	\$455
355 nm	355.0 nm	355.0 nm	364.9 – 900.0 nm	BLP01-355R-25	\$455
363.8 nm	363.8 nm	363.8 nm	374.0 – 900.0 nm	BLP01-364R-25	\$455
405 nm	400.0 nm	410.0 nm	421.5 – 900.0 nm	BLP01-405R-25 ★	\$405
441.6 nm	441.6 nm	441.6 nm	454.0 – 900.0 nm	BLP01-442R-25	\$405
457.9 nm	439.0 nm	457.9 nm	470.7 – 900.0 nm	BLP01-458R-25	\$405
473 nm	473.0 nm	473.0 nm	486.2 – 900.0 nm	BLP01-473R-25	\$405
488 nm	486.0 nm	491.0 nm	504.7 – 900.0 nm	BLP01-488R-25 ★	\$405
514.5 nm	505.0 nm	515.0 nm	529.4 – 900.0 nm	BLP01-514R-25	\$405
532 nm	532.0 nm	532.0 nm	546.9 – 900.0 nm	BLP01-532R-25 ★	\$405
561.4 nm	561.4 nm	561.4 nm	577.1 – 900.0 nm	BLP02-561R-25	\$405
568.2 nm	561.4 nm	568.2 nm	584.1 – 900.0 nm	BLP01-568R-25	\$405
594 nm	593.5 nm	594.3 nm	610.9 – 900.0 nm	BLP01-594R-25	\$405
632.8 nm	632.8 nm	632.8 nm	650.5 – 1200.0 nm	BLP01-633R-25 ★	\$405
635 nm	632.8 nm	642.0 nm	660.0 – 1200.0 nm	BLP01-635R-25	\$405
647.1 nm	647.1 nm	647.1 nm	665.2 – 1200.0 nm	BLP01-647R-25	\$405
664 nm	664.0 nm	664.0 nm	682.6 – 1200.0 nm	BLP01-664R-25	\$405
785 nm	780.0 nm	790.0 nm	812.1 – 1200.0 nm	BLP01-785R-25 ★	\$405
808 nm	808.0 nm	808.0 nm	830.6 – 1600.0 nm	BLP01-808R-25	\$405
830 nm	830.0 nm	830.0 nm	853.2 – 1600.0 nm	BLP01-830R-25	\$405
980 nm	980.0 nm	980.0 nm	1007.4 – 1600.0 nm	BLP01-980R-25	\$455
1064 nm	1064.0 nm	1064.0 nm	1093.8 – 1600.0 nm	BLP01-1064R-25 ★	\$455
1319 nm	1319.0 nm	1319.0 nm	1355.9 – 2000.0 nm	BLP02-1319R-25	\$455
1550 nm	1550.0 nm	1550.0 nm	1593.4 – 2000.0 nm	BLP01-1550R-25 ★	\$455

## Short Wave Pass

Nominal Laser Wavelength	Laser Wavelength Range		Passband	Part Number	Price
	$\lambda$ short	$\lambda$ long			
532 nm	532.0 nm	532.0 nm	350.0 – 517.1 nm	BSP01-532R-25	\$455
632.8 nm	632.8 nm	647.1 nm	350.0 – 615.1 nm	BSP01-633R-25	\$455
785 nm	780.0 nm	790.0 nm	350.0 – 758.2 nm	BSP01-785R-25	\$455

★ = High Selling Item – quantity ≤ 3 usually ships in 3 days, demand dependent.

 SEE SPECTRAL GRAPHS AND GET NUMERIC DATA FOR ALL OUR FILTERS AT [IDEX-HS.COM/SEMROCK](http://IDEX-HS.COM/SEMROCK)

# Edge Filters vs. Notch Filters for Raman Instrumentation

## RazorEdge® Filter Advantages:




- › Steepest possible edge for looking at the smallest Stokes shifts
- › Largest blocking of the laser line for maximum laser rejection

## StopLine® Notch Filter Advantages:

- › Measure Stokes and anti-Stokes signals simultaneously
- › Greater angle-tunability and bandwidth for use with variable laser lines

The graph on the top illustrates the ability of a long-wave-pass (LWP) filter to get extremely close to the laser line. The graph in the center compares the steepness of an edge filter to that of a notch filter. A steeper edge means a narrower transition width from the laser line to the high-transmission region of the filter. With transition widths guaranteed to be below 1% of the laser wavelength (on Semrock U-grade edge filters), these filters don't need to be angle-tuned!

The graph on the bottom shows the relative tuning ranges that can be achieved for edge filters and notch filters. Semrock edge filters can be tuned up to 0.3% of the laser wavelength. The filters shift toward shorter wavelengths as the angle of incidence is increased from 0° to about 8°. Semrock notch filters can be tuned up to 1.0% of the laser wavelength. These filters also shift toward shorter wavelengths as the angle of incidence is increased from 0° up to about 14°.

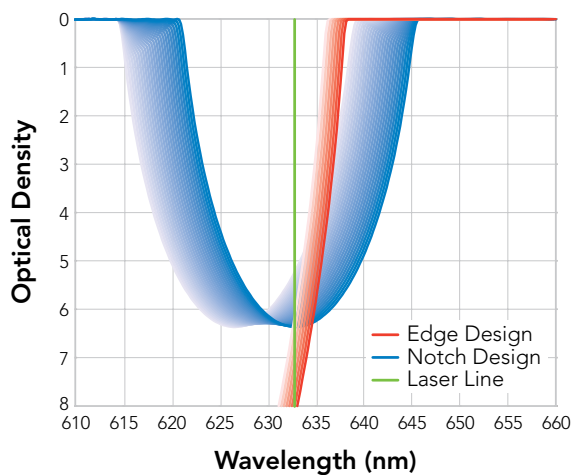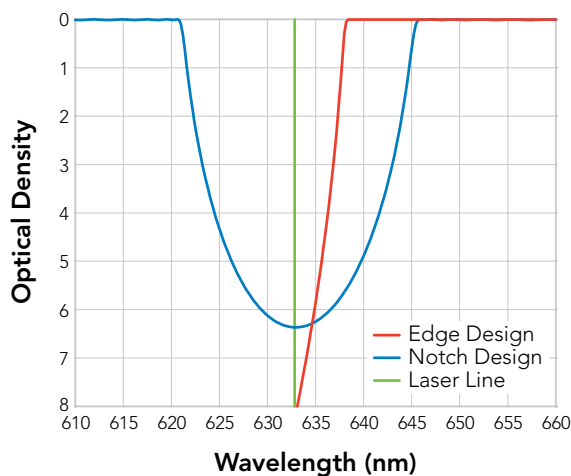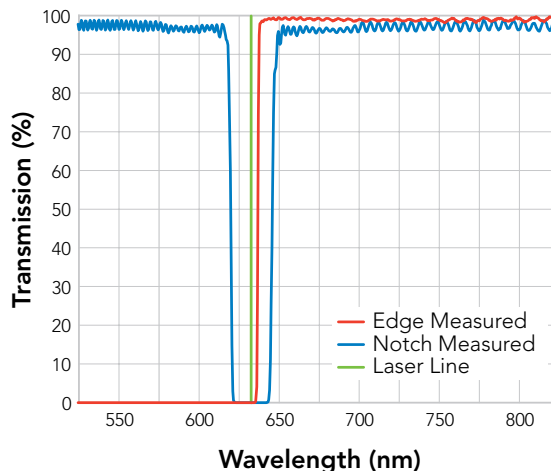




# StopLine® Single-notch Filters

StopLine deep notch filters rival the performance of holographic notch filters but in a less expensive, more convenient, and more reliable thin-film filter format (U.S. Patents No. 7,123,416 and No. 9,354,370). These filters are ideal for applications including Raman spectroscopy, laser-based fluorescence instruments, and biomedical laser systems.


- › The stunning StopLine E-grade notch filters offer high transmission over ultra-wide passbands (UV to 1600 nm)
- › Deep laser-line blocking for maximum laser rejection (OD > 6)
- › High laser damage threshold and proven reliability
- › Rejected light is reflected for convenient alignment and best stray-light control
- › Multi-notch filters are available for blocking multiple laser lines (visit [idex-hs.com](http://idex-hs.com) for details)



Semrock introduced a breakthrough invention in thin-film optical filters: our StopLine E-grade thin-film notch filters have ultrawide passbands with deep and narrow laser-line blocking. Unheard of previously in a thin-film notch filter made with multiple, discrete layers, these notch filters attenuate the laser wavelength with OD > 6 while passing light from the UV well into the near-infrared (1600 nm). They are especially suited for optical systems addressing multiple regions of the optical spectrum (e.g., UV, Visible, and Near-IR), and for systems based on multiple detection modes (e.g., fluorescence, Raman spectroscopy, laser-induced breakdown spectroscopy, etc.).

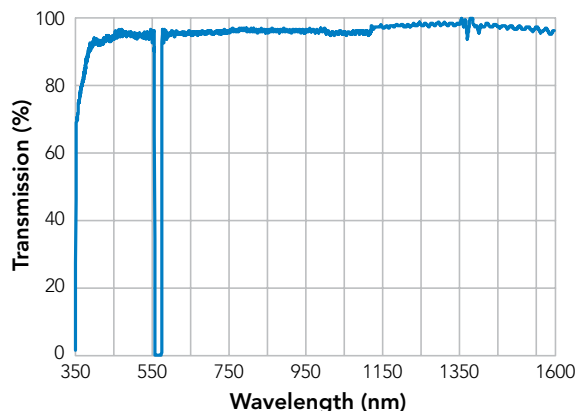

Typical measured data for NF03-561E.

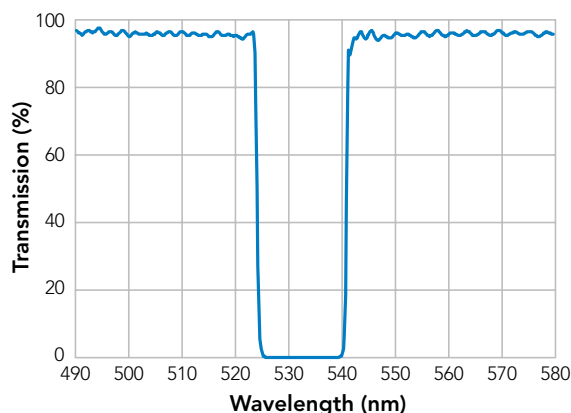


Typical measured data for NF03-532E.

## StopLine Single-notch Filters

Wavelength	Passband Range	Typical 50% Notch Bandwidth	Laser-line Blocking	Part Number	Price
405.0 nm	330.0 – 1600.0 nm	9 nm	OD > 6	NF03-405E-25	\$945
473.0 nm	350.0 - 1600.0 nm	13 nm	OD > 6	NF03-473E-25	\$945
488.0 nm	350.0 – 1600.0 nm	14 nm	OD > 6	NF03-488E-25 ★	\$945
514.5 nm	350.0 – 1600.0 nm	16 nm	OD > 6	NF03-514E-25	\$990
532.0 nm	350.0 – 1600.0 nm	17 nm	OD > 6	NF03-532E-25	\$945
532.0 nm	399.0 – 709.3 nm	17 nm	OD > 6	NF01-532U-25 ★	\$765
561.4 nm	350.0 – 1600.0 nm	19 nm	OD > 6	NF03-561E-25	\$990
594.1 nm	350.0 – 1600.0 nm	22 nm	OD > 6	NF03-594E-25	\$945
632.8 nm	350.0 – 1600.0 nm	25 nm	OD > 6	NF03-633E-25 ★	\$945
642.0 nm	350.0 - 1600.0 nm	26 nm	OD > 6	NF03-642E-25	\$990
658.0 nm	350.0 – 1600.0 nm	27 nm	OD > 6	NF03-658E-25	\$990
785.0 nm	350.0 – 1600.0 nm	39 nm	OD > 6	NF03-785E-25 ★	\$945
808.0 nm	350.0 – 1600.0 nm	41 nm	OD > 6	NF03-808E-25	\$945

★ = High Selling Item – quantity ≤ 3 usually ships in 3 days, demand dependent.

 **LOOKING FOR A 1064 nm NOTCH FILTER?**  
TRY THE NF03-532/1064E AT [IDEX-HS.COM](http://IDEX-HS.COM)

## Q MEASURE DEEPER BLOCKING

# KolaDeep™ Spectral Measurement System

New fluorescence-based life science and biomedical instrumentation increasingly requires high-performance optical filters with very high blocking (OD) and steep spectral edges that transition between high transmission and high OD. Our IDEX Health & Science Semrock engineers therefore developed the KolaDeep Spectral Measurement System (SMS), a proprietary award-winning platform for proven measurements of the steepest and deepest spectral features of our optical filters. This system is used in routine production to ensure that instruments made from these filters will deliver their intended performance.


The platform's capabilities include:

- ▶ KolaDeep measures blocking to OD 8-9 across wavelength spans in the UV, visible and NIR ranges, and spectral features to OD 11
- ▶ KolaDeep resolves edges steeper than 0.2% relative to the edge wavelength from 90% transmission to OD > 7


Our engineering team performed extensive qualification testing on KolaDeep to characterize and confirm the wavelength and photometric accuracy of the system

The recent white paper "Advanced Spectral Measurement Systems at IDEX Health & Science Semrock" presents examples and analyses of Semrock filters based on measurements by the KolaDeep SMS; a few are shown here


The graph below shows OD measured on a custom penta-band filter. The blocking between passbands is specified to exceed OD 8 to prevent excitation light leakage. Theory and measurement agree to OD > 8, with some blocking reaching from OD > 9 to the noise floor.

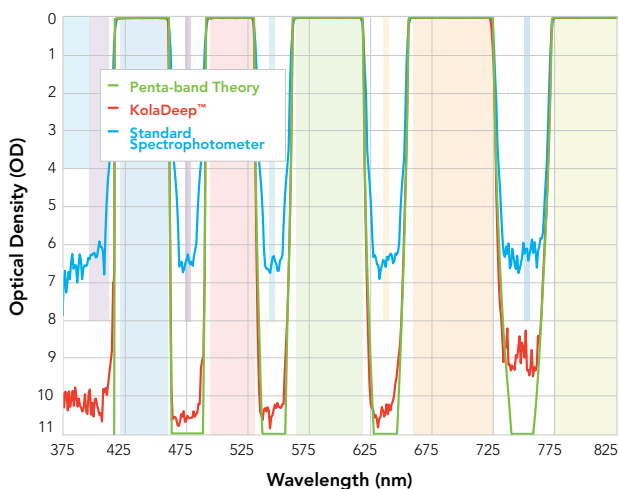

The graph below shows spectral performance for a custom hexa-notch filter designed to ensure precise blocking of six laser lines. This customer required assurance that each laser line is repeatably blocked at a specific level. The spectra measured with our KolaDeep SMS show precise matches to the specifications, even in the UV region. The measured blocking (red) matches theory (green).

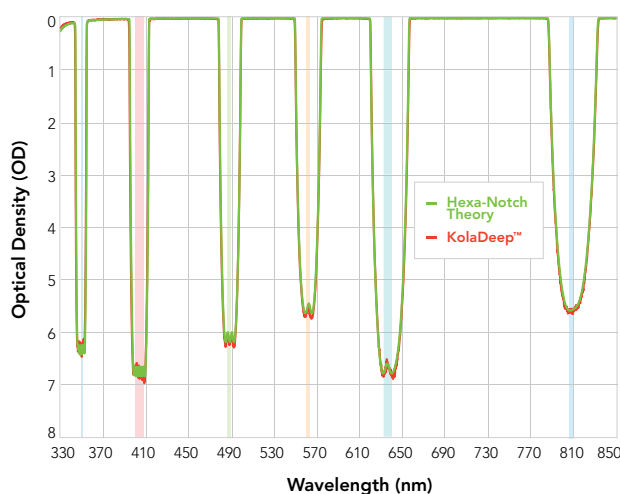

The KolaDeep SMS performs well over the UV and NIR spectral range. An example is shown below for the 1064 nm Semrock MaxLine® laser-line cleanup filter. Measurements (red) confirm the required blocking and agree with the theoretical design (green) to OD > 7. Data (blue) from a superior commercially available system are not accurate beyond OD 3.

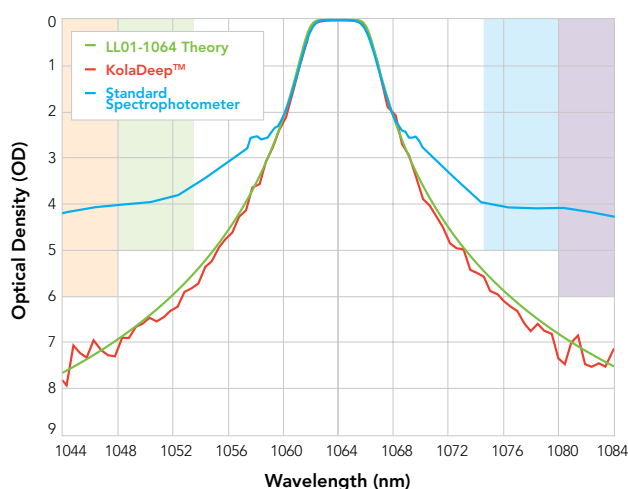

## TECHNICAL NOTE


# Ultraviolet (UV) Raman Spectroscopy

Raman spectroscopy measurements generally face two limitations: (1) Raman scattering cross sections are tiny, requiring intense lasers and sensitive detection systems just to achieve enough signal; and (2) the signal-to-noise ratio is further limited by fundamental, intrinsic noise sources like sample autofluorescence. Raman measurements are most commonly performed with green, red, or near-infrared (IR) lasers, largely because of the availability of established lasers and detectors at these wavelengths. However, by measuring Raman spectra in the ultraviolet (UV) wavelength range, both of the above limitations can be substantially alleviated.

Visible and near-IR lasers have photon energies below the first electronic transitions of most molecules. However, when the photon energy of the laser lies within the electronic spectrum of a molecule, as is the case for UV lasers and most molecules, the intensity of Raman-active vibrations can increase by many orders of magnitude – this effect is called “resonance-enhanced Raman scattering.”

Although UV lasers tend to excite strong autofluorescence, it typically occurs only at wavelengths above about 300 nm, independent of the UV laser wavelength. Since even a  $4000\text{ cm}^{-1}$  (very large) Stokes shift leads to Raman emission below 300 nm when excited by a common 266 nm laser, autofluorescence simply does not interfere with the Raman signal making high signal-to-noise ratio measurements possible.

An increasing number of compact, affordable, and high-power UV lasers have become widely available, such as quadrupled, diode-pumped Nd:YAG lasers at 266 nm and NeCu hollow-cathode metal-ion lasers at 248.6 nm, making ultra-sensitive UV Raman spectroscopy a now widely accessible technique.

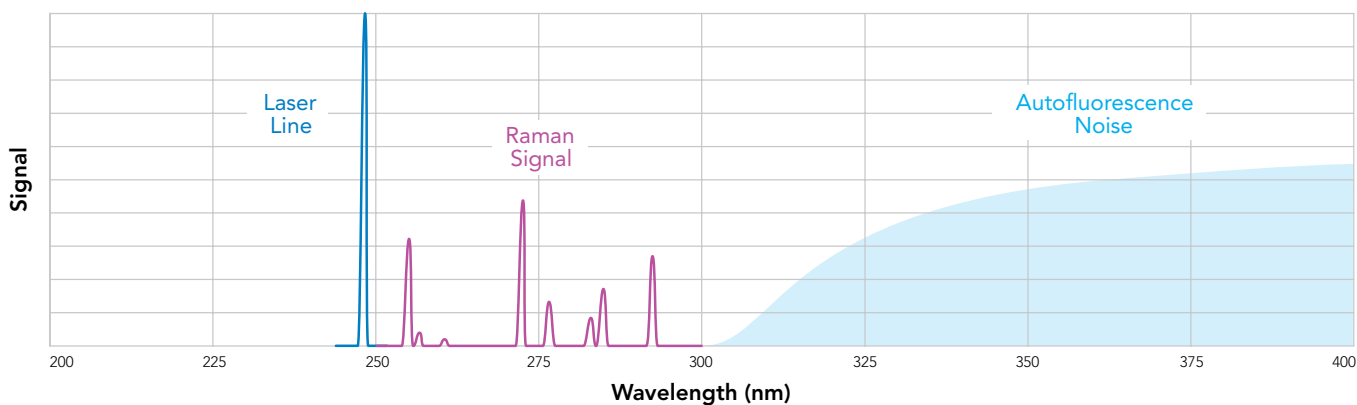


## TECHNICAL NOTE

# Coherent Raman Scattering (CRS, CARS and SRS)

With coherent Raman scattering (CRS) it is possible to perform highly specific, label-free chemical and biological imaging with orders of magnitude higher sensitivity at video-rate speeds compared with traditional Raman imaging. CRS is a nonlinear four-wave mixing process that is used to enhance the weak spontaneous Raman signal associated with specific molecular vibrations. Two different types of CRS that are exploited for chemical and biological imaging are coherent anti-Stokes Raman scattering (CARS) and stimulated Raman scattering (SRS).

In CRS, two lasers are used to excite the sample. The wavelength of a first laser (often a fixed-wavelength, 1064 nm laser) is set at the Stokes frequency,  $\omega_{\text{Stokes}}$ . The wavelength of the second laser is tuned to the pump frequency,  $\omega_{\text{pump}}$ . When the frequency difference  $\omega_{\text{pump}} - \omega_{\text{Stokes}}$  between these two lasers matches an intrinsic molecular vibration of frequency  $\Omega$  both CARS and SRS signals are generated within the sample.


In CARS, the coherent Raman signal is generated at a new, third wavelength, given by the anti-Stokes frequency  $\omega_{\text{CARS}} = 2\omega_{\text{pump}} - \omega_{\text{Stokes}} = \omega_{\text{pump}} + \Omega$ . In SRS there is no signal at a wavelength that is different from the laser excitation wavelengths. Instead, the intensity of the scattered light at the pump wavelength experiences a stimulated Raman loss (SRL), with the intensity of the scattered light at the Stokes wavelength experiencing a stimulated Raman gain (SRG). The key advantage of SRS microscopy over CARS microscopy is that it provides background-free chemical imaging with improved image contrast, both of which are important for biomedical imaging applications where water represents the predominant source of nonresonant background signal in the sample.

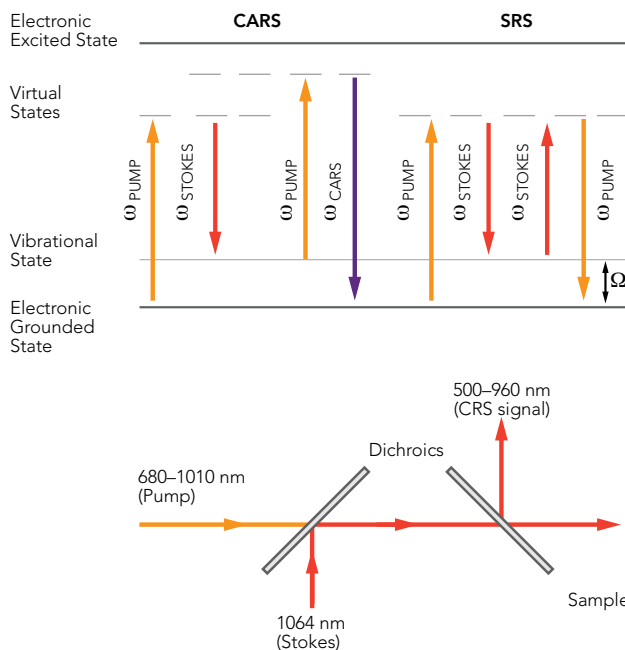

Coherent Raman scattering energy diagrams for both CARS and SRS (top), and a schematic of a typical experimental setup (bottom).

## Harmonic Generation Microscopy

Harmonic generation microscopy (HGM) is a label-free imaging technique that uses high-peak power ultrafast lasers to generate appreciable image contrast in biological imaging applications. Harmonic generation microscopy exploits intrinsic energy-conserving second and third order nonlinear optical effects. In second-harmonic generation (SHG) two incident photons interact at the sample to create a single emission photon having twice the energy i.e.,  $2\omega_1 = \omega_{\text{SHG}}$ . A prerequisite for SHG microscopy is that the sample must exhibit a significant degree of noncentrosymmetric order at the molecular level before an appreciable SHG signal can be generated. In third-harmonic generation (THG), three incident photons interact at the sample to create a single emission photon having three times the energy i.e.,  $3\omega_1 = \omega_{\text{THG}}$ . Both SHG and THG imaging techniques can be combined with other nonlinear optical imaging (NLO) modalities, such as multiphoton fluorescence and coherent Raman scattering imaging. Such a multimodal approach to biological imaging allows a comprehensive analysis of a wide variety of biological entities, such as individual cells, lipids, collagen fibrils, and the integrity of cell membranes at the same time.

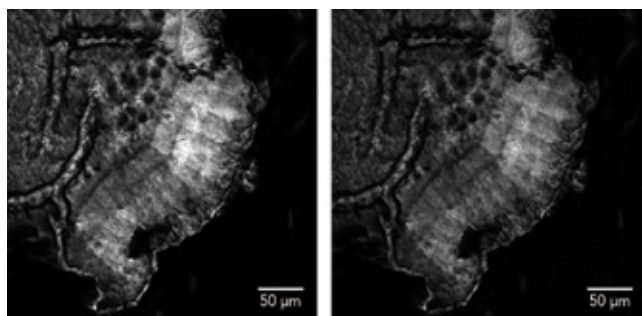

Coherent anti-Stokes Raman (CARS) imaging of cholesteryl palmitate. The image on the left was obtained using Semrock filter FF01-625/90. The image on the right was obtained using a fluorescence bandpass filter having a center wavelength of 650 nm and extended blocking. An analysis of the images revealed that the FF01-625/90 filter provided greater than 2.6 times CARS signal. Images courtesy of Prof. Eric Potma (UC Irvine).

## TECHNICAL NOTE

# Filter Spectra at Non-normal Angles of Incidence

While most applications call for optical filters to be used at normal incidence, it is important to understand how the spectral properties of different types of filters change when using these filters at non-normal angles of incidence (AOI). There are two main effects exhibited by all filter spectra as the angle is increased from normal:

1. The features of the spectrum shift to shorter wavelengths.
2. Two distinct spectra emerge – one for s-polarized light and one for p-polarized light.

## RazorEdge® Edge Filters

The graph in Figure 1 shows a series of spectra derived from a typical RazorEdge long-wave-pass (LWP) filter design, and can be applied to any of the RazorEdge edge filters. Here, the resulting wavelength  $\lambda$  is compared to the wavelength  $\lambda_0$  of the edge location at normal incidence. As the angle increases from normal incidence, the filter edge shifts toward shorter wavelengths and the edges associated with s- and p-polarized light shift by different amounts. For LWP filters, the response of p-polarized light shifts more than s-polarized light. The opposite is true for SWP filters. This polarization splitting causes the unpolarized spectrum to show a "shelf" near the 50% transmission point, but note that the edge steepness remains intact, even for polarized light.

The shift of any spectral feature of a Semrock filter can be calculated using the MyLight tool on the website-based product page. For example, on the product for the 785 nm RazorEdge long pass filter at 785 nm, click on the MyLight link in the [green box](#). Visit [idex-hs.com/resources/tools-drawings/mylight-tutorial](https://www.idex-hs.com/resources/tools-drawings/mylight-tutorial) for a tutorial for how to use MyLight. Also note that any catalog filter can also be modeled at specific AOI, CHA, and Polarization at [idex-hs.com/SearchLight](https://www.idex-hs.com/SearchLight).

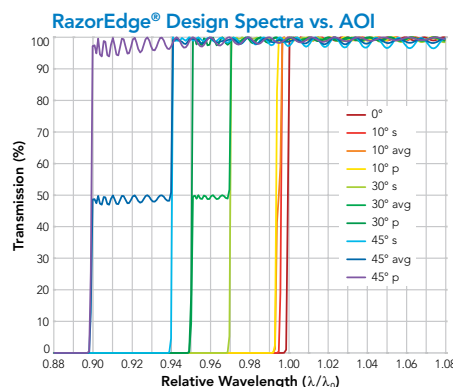


Figure 1: Two distinct filter spectra emerge for s- and p-polarized light as AOI increases.

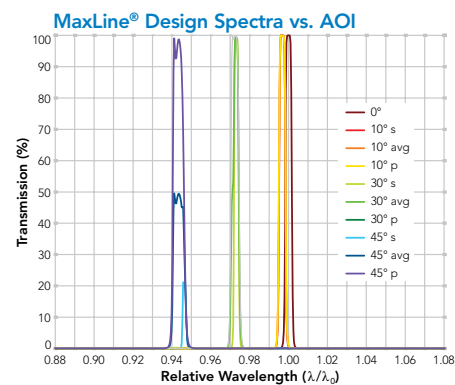


Figure 2: MaxLine filter spectral response to AOI shift.

## MaxLine® Laser-line Filters

Varying filter designs will respond to changes in AOI in different ways, but still show marked differences between s- and p-polarized light. The spectra of a MaxLine laser-line filter (Figure 2) shows that as the angle increases from normal incidence, the center wavelength shifts toward shorter wavelengths and the bandwidth broadens slightly for p-polarized light while narrowing for s-polarized light. The most striking feature is the decrease in transmission for s-polarized

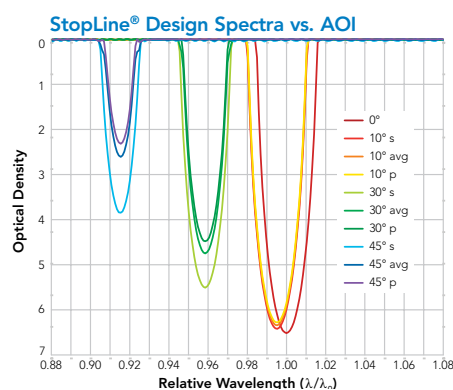


Figure 3: U- and S-grade StopLine filter spectra.

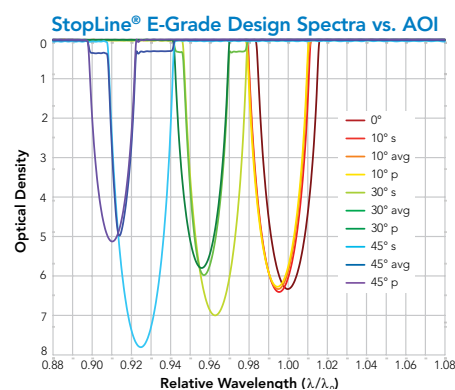


Figure 4: E-grade StopLine filter spectra.

light, while the p-polarized light maintains high transmission. Use the MyLight Tool described above to model and download spectra of AOI-shifted MaxLine filters, or model them in [SearchLight](#).

## StopLine® Notch Filters

The StopLine series of filters has a different response to AOI depending on the grade of filter used. For the standard U- and S-grade filters (Figure 3), as the angle is increased from normal incidence, the notch center wavelength shifts to shorter wavelengths, the notch depth decreases, and the notch bandwidth decreases (with a greater decrease for p-polarized light than for s-polarized light).

The E-grade filter (Figure 4), which is based on a different type of filter design, shows a large increase in OD for s-polarized light as the AOI increases. As the angle is increased from normal incidence, the notch center wavelength shifts to shorter wavelengths; however, the shift is greater for p-polarized light than it is for s-polarized light.



To learn more about our Semrock notch filters, search for StopLine® Notch Filters at [idex-hs.com](https://www.idex-hs.com). Use our MyLight Tool to model and download spectra of AOI-shifted StopLine filters, or model them in [SearchLight](#).



Semrock, Inc. manufactures optical filters that set the standard for the life science and analytical instrumentation industries, as well as optical filters and mirrors for laser and optical systems. OEM filters are manufactured in volume using thin-film sputtering and proprietary volume manufacturing technology. All Semrock products carry a ten-year warranty.

Founded in September 2000, Semrock is based in Rochester, New York, a well-known world center for optics, and has sales offices throughout the United States. In October 2008 Semrock became a unit of IDEX Corporation.

### Address:

IDEX Health & Science, LLC  
Center of Excellence  
1180 John Street  
Rochester, NY 14586  
USA

### OEM Customers North America:

To submit purchase orders, email:  
IHSOrders@idexcorp.com

For technical questions, email:  
semrock@idexcorp.com

### Non-OEM Customers North America:

Visit:  
www.avr-optics.com/semrock

Email:  
info@avr-optics.com

### Online Ordering:

www.idex-hs.com

### Phone:

Worldwide:  
+1 585 625 5000

Within US  
and Canada:  
866 736 7625  
(866-SEMROCK)

### General Inquiries:

Phone:  
+1 585 594 7000

Email:  
semrock@idexcorp.com

Fax:  
+1 585 594 7095

## Ordering

- **OEM customers in the USA and Canada** may place orders on the IDEX Health & Science website using a credit card. An account is required to order — visit [www.idex-hs.com](http://www.idex-hs.com) for details on setting up your account. Please note that credit approval is subject to review.
- **North American Government, University, and National Laboratory customers** may purchase Semrock Optical Filters through our domestic distribution partner, AVR Optics (account required, [www.avr-optics.com](http://www.avr-optics.com)).
- **International customers** (outside the USA and Canada) may purchase Semrock Optical Filters through Semrock's distribution network. Please visit our website, [www.idex-hs.com](http://www.idex-hs.com), to find your local distributor.

## Online Quotes

- OEM Customers in the USA and Canada can generate quotes online (account required, [www.idex-hs.com](http://www.idex-hs.com)).
- North American Government, University, and National Laboratory customers can generate quotations through our domestic distribution partner, AVR Optics (account required, [www.avr-optics.com](http://www.avr-optics.com)).
- International customers purchasing through Semrock's distribution network may send shopping cart contents to their local distributor (account required, [www.idex-hs.com](http://www.idex-hs.com)).

## Pricing & Availability

All prices are domestic USD and subject to change without notice. Check our website for current pricing and availability.

## Shipping

Orders received for catalog products will be reviewed for best possible delivery date within 48 hours. Domestic orders are shipped UPS Ground or 2-day, unless otherwise requested. We also accept customers' shipping account numbers for direct carrier billing.

## 10-Year Warranty

Be confident in your filter purchase with our comprehensive 10-year warranty. Built to preserve their high level of performance in test after test, year after year, our filters reduce your cost of ownership by eliminating the expense and uncertainty of replacement costs.

## 30-Day Return Policy

Semrock offers a 30-day return policy for unused, undamaged standard-sized catalog products in their original packaging. Returns must be made within 30 days of purchase with an RMA number to qualify for a full credit. OEM volumes returned after 30 days may incur a restocking fee, and custom-sized parts are not eligible for return. Products purchased through distributors must be returned via the distributor, and all returns are subject to inspection. Returns without an RMA number will not be accepted and will be returned at the sender's expense.

## RoHS & REACH Compliant

Semrock complies with both RoHS and REACH regulations. You can find our full compliance statement at [www.idex-hs.com/about/environmental-policy](http://www.idex-hs.com/about/environmental-policy)

## ISO 9001:2015 Certified

Semrock's Quality management system is certified to ISO 9001:2015.



ISO 9001:2015

www.tuv.com  
ID 9108657280